

 HENRY SCHEIN®

ORTHODONTICS

Product Exclusive




smartIPRTM

The Smarter Way to Perform IPR

smartIPR 


Safer IPR Starts Here

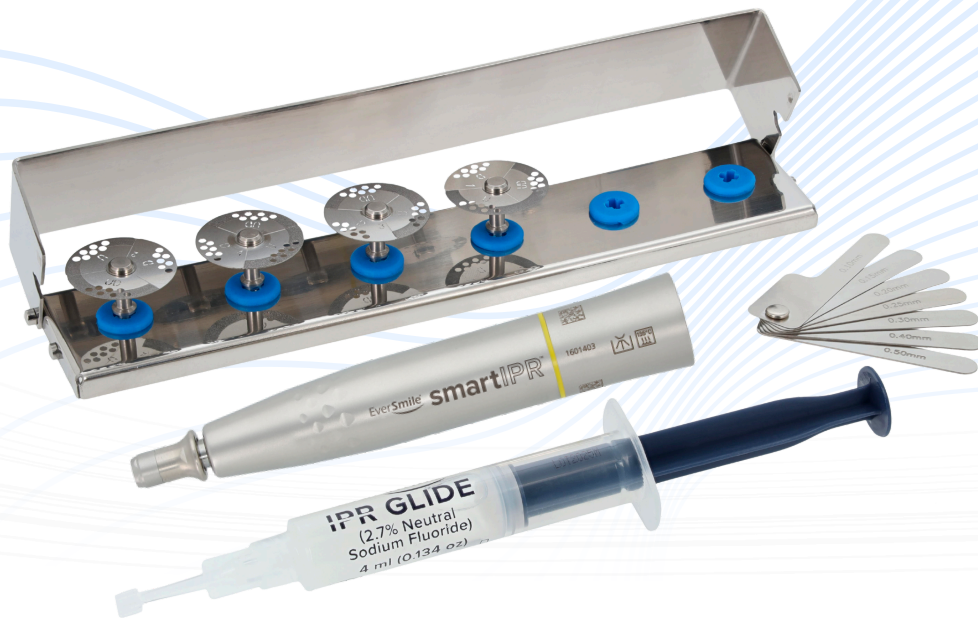

The **Eversmile® smartIPR™ system** advantage is its ability to open tight interproximal contacts gradually with precision, eliminating the need for traditional, razor-sharp hand strips or spinning disks.

The Benefits of smartIPR

The smartIPR system delivers safe, precise space creation by polishing enamel surfaces – reducing the risk of cutting hard and soft tissue.

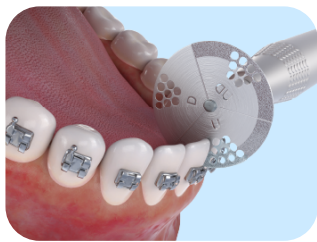


Precision Without Rotation

Combining innovative design with **oscillation, vibration, and lubrication**, smartIPR redefines IPR.

Advanced Enamel Protection

Eversmile's IPRGlide™ enhances IPR procedures with **lubrication and Rx fluoride** delivery to abraded enamel surfaces.

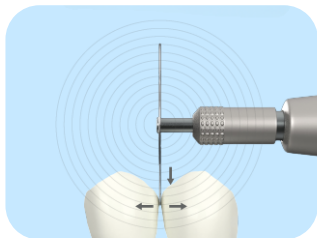

Protection by Design


smartIPR Disks create space through **vibrating, gentle oscillation** while helping polish enamel surfaces.

[Click here](#)
to see smartIPR
in action.



Why smartIPR is the Smart Choice

Compare the smartIPR System to traditional IPR strips and discs to see how precision-guided technology can improve safety, efficiency, and workflow control.

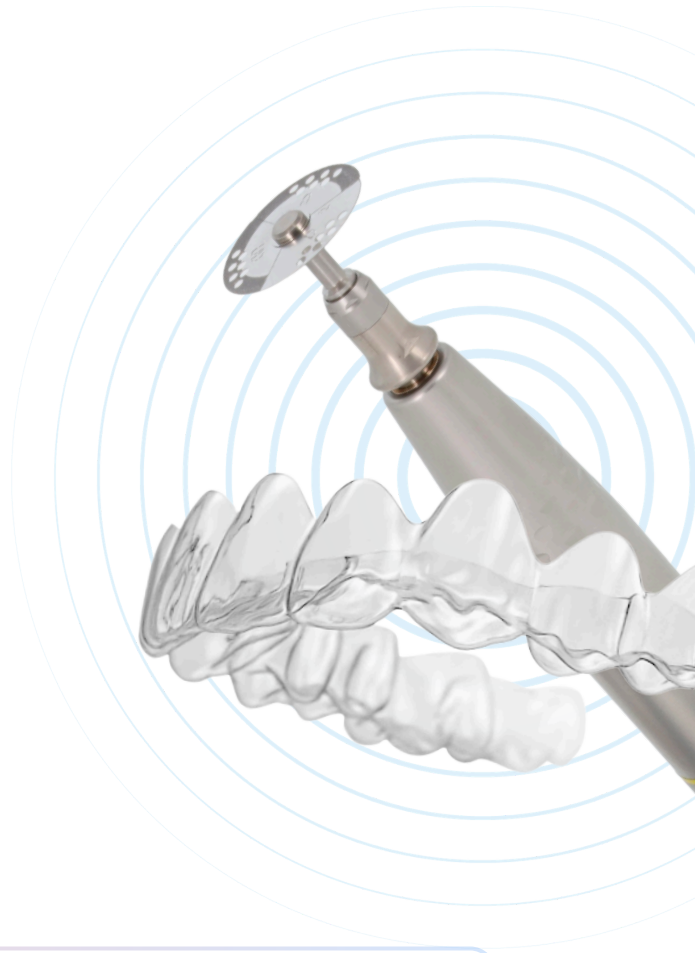

Traditional IPR Strips

- Manual hand stripping
- Variable reduction
- Soft tissue contact risk
- Multiple strip sizes


Traditional IPR Discs

- Freehand disc reduction
- Variable reduction
- Adjacent tissue risk
- Multiple discs & mandrels

smartIPR SYSTEM

- Precision-guided workflow
- Controlled reduction
- Designed for soft tissue safety
- 3-in-1 disk system

The smartIPR System Includes

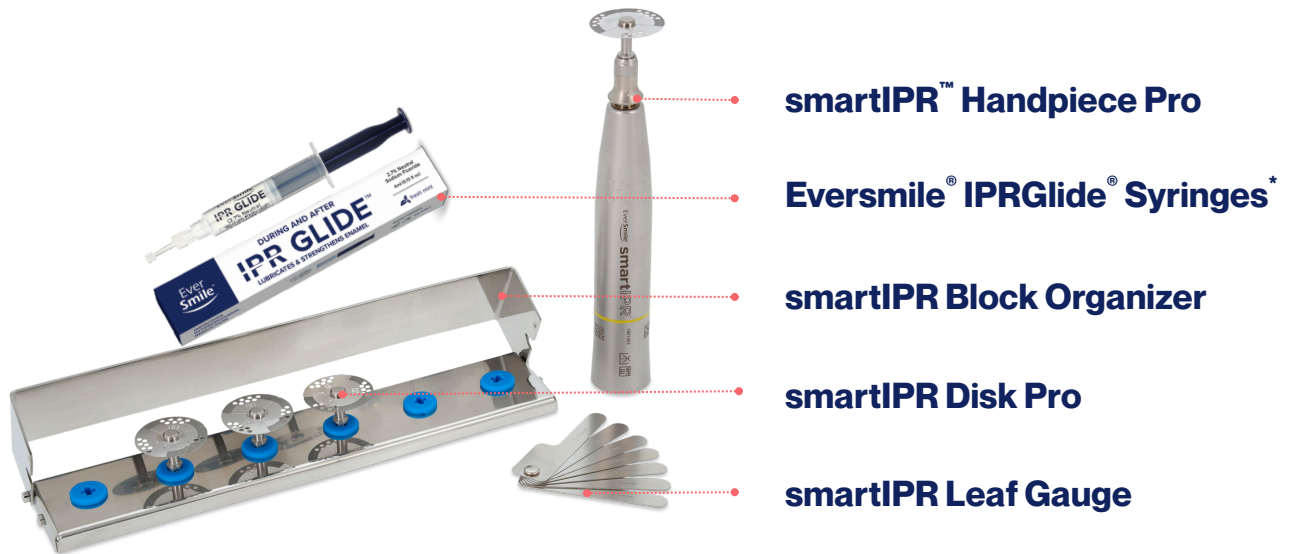

Compatible Motor Types:

- **Electric:** Many electric low-speed motors use short E-type connection. Recommended: 20,000–30,000 RPM
- **Air-Driven:** Use air-drive motors with a short E-type connection. Recommended: Up to 20,000 RPM

[Click here](#) to view the smartIPR Quickstart Guide.



Coming Soon: smartIPR Air Motor, Short E-Type**

Check with your sales representative for availability.




*Eversmile IPRGlide Syringes are not included in the smartIPR Kit Pro.


**smartIPR Air Motor coming soon. Check with your sales representative for availability.

Smarter 3-in-1 Disk Design

The smartIPR Disk Pro features a patented non-spinning design with three working sectors built into each disk, providing superior control and safeguard adjacent tooth surfaces, ensuring reduction occurs only where intended.

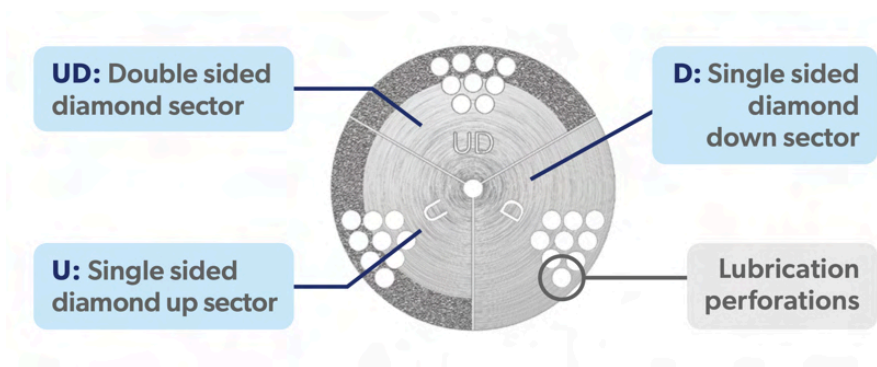

4 Disk Sizes + 3 Sectors Each = 12 Instruments

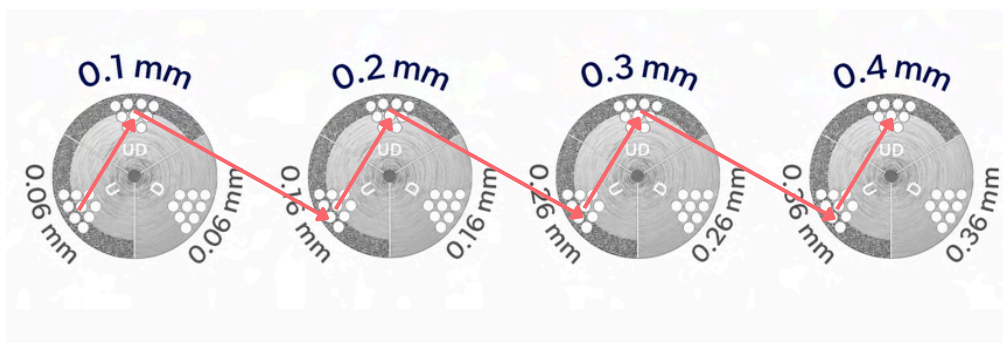

Pro Tip: The 0.1 mm and 0.2 mm smartIPR Disks are commonly used; maintain backup inventory and replace worn disks as needed.

Controlled Vibration. Enhanced Protection.




Precision Without Rotation

The smartIPR Handpiece Pro delivers precise, controlled vibration instead of spinning motion. This innovative design helps protect hard and soft tissue during treatment.

Sonic Vibration at Work

Powered by sonic vibration, smartIPR Disks operate at 128–192 Hz — up to 100 times faster than traditional hand-strip motion.*

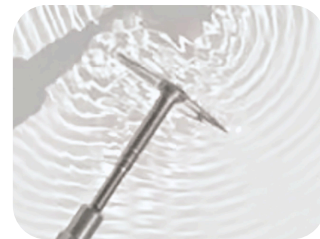

The Smarter Way to Create Space

By transferring energy through controlled vibration to the periodontal ligaments, smartIPR Disks gently move through contacts to create space and polish enamel without cutting.

[Click here](#)
to see smartIPR
in action.



*Dependent on slow-speed motor settings

The smartIPR System

 **HENRY SCHEIN**
ORTHODONTICS

How to Use smartIPR

A key advantage of the smartIPR system is its ability to open tight interproximal contacts gradually and with precision, eliminating the need for traditional, razor-sharp hand strips.

1

Begin with the thinnest disk (0.1 mm) using one of the single-sided segments ("U" or "D"), which provides an approximate working thickness of 0.065 mm.


2

Place IPRGlide on the smartIPR Disk or on the occlusal surface between the teeth you wish to perform IPR on.




3

Position the smartIPR Disk parallel to the teeth's roots.

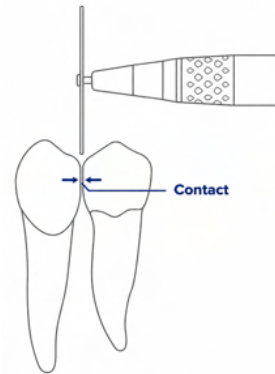

4

Oscillate the lubricated disk and advance it apically into the interproximal contact. The oscillatory motion transfers vibrational energy to the teeth and periodontal ligaments, facilitating slight separation. This permits the disk to pass through the contact area, while gently removing enamel, and creating the desired opening.

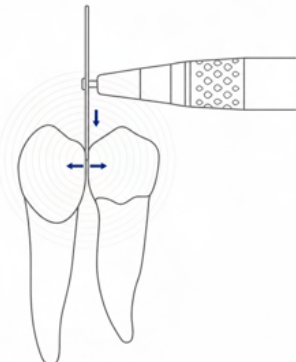


[Click here](#) to view the smartIPR Quickstart Guide.

How to Use smartIPR

5

Apply light pressure to move the oscillating disk apically with a gentle painting motion. Sweep the disk buccal-lingually and apico-incisally until the smartIPR Disk passes through the contact. Continue refining the surface(s) until the next size smartIPR Disk fits into the space.

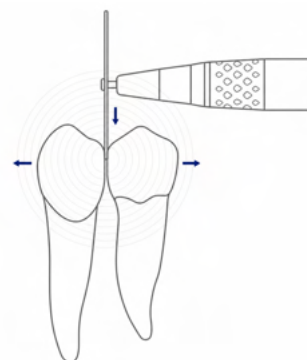

6

After the contact is gently opened, rotate to the double-coated segment ("UD") to increase the space to the full 0.1 mm.

7

Once the initial contact has been opened with the 0.1 mm smartIPR Disk, continue sequentially with the 0.2 mm disk. Start again with the "U" or "D"-single-coated surface (approximately 0.165 mm), then transition to the "UD" double-coated 0.2 mm segment.


8

Continue progressing in this stepwise manner to achieve the desired amount of interproximal reduction.

[Click here](#) to view
the smartIPR
Quickstart Guide.



The smartIPR Ordering Guide

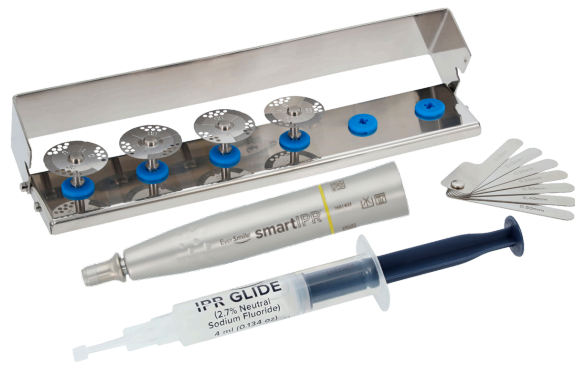

DESCRIPTION	QTY	HSD ITEM #
smartIPR Kit Pro*  Includes: smartIPR Handpiece Pro (1), 4 smartIPR Disk Pro (1 of each 0.1mm, 0.2mm, 0.3mm, 0.4mm), smartIPR Disk Block (1), smartIPR Leaf Gauge (1)	1 KIT	150-9227
smartIPR Handpiece Pro	1	151-1498
smartIPR Disk Pro 0.1 mm	1	151-1492
smartIPR Disk Pro 0.2 mm	1	151-1493
smartIPR Disk Pro 0.3 mm	1	151-1494
smartIPR Disk Pro 0.4 mm	1	151-1495
smartIPR Disk Block	1	151-1497
smartIPR Leaf Gauge	1	151-1496
smartIPR Air Motor, Short E-Type**	1	Coming Soon
Eversmile IPRGlide Syringe**	1	Coming Soon



Pro Tip: Since they are commonly used, we recommend purchasing an additional 3 qty 0.1 mm and 1 qty 0.2 mm smartIPR Disks with your Kit Pro purchase.

*Eversmile IPRGlide Syringes are not included in the smartIPR Kit Pro.

**smartIPR Air Motor and Eversmile IPRGlide Syringes coming soon.
smartIPR system is not available in all countries. See your local distributor for details.

Your smartIPR Resources

Visit the smartIPR Resource Hub for clinical guides, workflow videos, and more!



View helpful online smartIPR guides & videos online



Watch smartIPR with real patient applications



Download the smartIPR Quick Start Guide



Review smartIPR chairside workflow guides & videos



Signup for upcoming clinical webinars

[Visit Resources.HenryScheinOrtho.com](https://www.henryscheinortho.com/smartIPR)



Start Loving IPR with the smartIPR System




FIND YOUR PERFECT ORTHODONTIC SOLUTION

Henry Schein® Orthodontics™ aims to deliver next generation orthodontic solutions and new treatment concepts to help you provide exceptional clinical outcomes, and to ensure that you, as an orthodontist, truly love what you do.



Exceed Your Expectations

Delivering new ideas and treatment concepts to help you provide exceptional clinical outcomes.

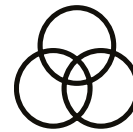

Beyond the Basics

We go beyond supplier-client relationship to cultivate partnerships with our customers to deliver positive patient outcomes.



State of the Art

Orthodontic products and solutions enable our customers to offer exceptional patient care.



World-Class Products

Products you have come to trust we create partnerships with you that go beyond the supplier-client relationship.



For more information on our products and educational offerings, please contact us:
HenryScheinOrtho.com
For phone orders, contact Customer Service: 888.851.0533 | 800.547.2000
To email an order: USASales@HenryScheinOrtho.com

Henry Schein Orthodontics, 800 King George Blvd, Suite 100 Savannah, GA 31419

© 2026 Henry Schein Orthodontics. All rights reserved. EverSmile, smartIPR System, EverSmile IPRGlide Syringes are registered trademarks of EverBrands, Inc. M3983 06/26